



NEW GFB LAB

A black and white photograph of a GFB LAB machine, a vertical industrial-scale laboratory mixer, is partially obscured by the large, bold, red and white text "NEW GFB LAB". The machine has a cylindrical body with various ports and a control panel on the right side.

GFB LAB

The New Benchmark for Product Development



EXPERIENCE THE NEW BENCHMARK FOR PRODUCT DEVELOPMENT





THE NEW GFB LAB POWDER DRYER-COATER-GRANULATOR FROM GLATT GSPL

*The modular multi-talent for full process
flexibility and product development*

The GFB LAB is the most versatile and adaptable system of all Glatt laboratory units. The flexible multi-function machine offers as many as three processing options in a single system: drying, granulating, coating. Apart from the drying or granulating insert, a Glatt HS Wurster insert for bottom spray are also available as options. From simple drying to demanding powder coating, no matter whether the top spray, bottom spray method is used – anything is possible.

The system is particularly easy to clean and is ideal for the production of clinical samples. A fluidized bed system operating in batch mode.

The GFB LAB is suitable for working volumes from 0.5 – 10 l.

Our Standard is the Result of our Claim:

FLEXIBILITY

- » Modular concept for maximum flexibility
- » Integrated deduster, de- humidification
- » Various process inserts for drying, granulating, Wurster coating

PRECISION

- » Increasing level of knowledge about processes (PAT)
- » Deeper process understanding (PAT)
- » Real time development (PAT)

CONTROL

- » Control system designed for easy and intuitive operation
- » Automated system in compliance with GAMP



OTHER HIGHLIGHTS

Designed for maximum safety

The closed system offers maximum pressure shock resistance for all types of dust (ST 1+2) and hybrid mixtures. Your operators, your products, and of course the environment are perfectly protected. Dust-ex and solvent-ex versions are available depending on your individual requirements.

The GFB LAB is equipped with an additional dehumidification module for dehumidification of inlet air down to a dew point of 10°C.

It fits smoothly into the design of the basic GFB LAB module. It is mobile on wheels and can be installed close to the machine or in the vicinity.



Ergonomic design for maximum efficiency


All functions are perfectly integrated into the design and enable optimal functionality: starting with flexible process inserts that can be pulled out and rotated, to removable product containers and through to well-thought-out connectors and a plug & play function to start-up right away. Easy charging and discharging as well as the user-friendly cleaning design further support efficient processing. The practical worktop made from synthetic granite is solvent, scratch, and impact resistant and will still look like new after many years of use.





Maximum process flexibility

The new GFB LAB offers the possibility to dry, granulate or coat particles in a single machine. This innovative concept allows to spray systems to be converted to bottom spray systems in just a few minutes. All with just one single product container. With the GFB LAB you are prepared for future requirements. Additionally needed process components can be retrofitted at any time at a low cost. This gives you a choice between three different process inserts for drying, granulating, and particle coating. The all-rounder available in two sizes 3 l and 6 l: Drying, granulating and coating with the unique Combo container. Top and bottom spray. The maximizer 12 l for LAB scale: Drying and top spray granulation with high volume product container

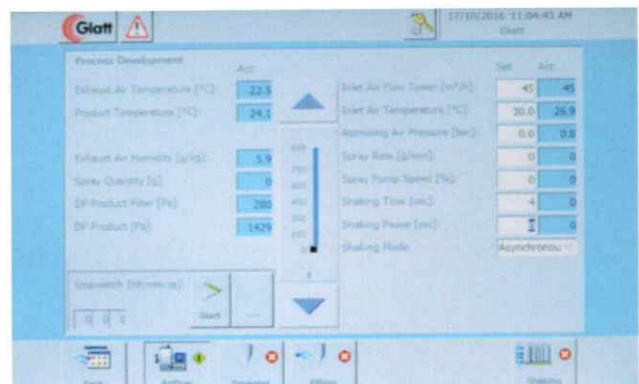

PAT for optimized process control

The integration of a Process Analytical Technology (PAT) system greatly enhances product quality by monitoring and controlling a whole suite of process parameters, such as particle size, moisture and blend uniformity. The smart design of the GFB Lab's product bowl in combination with the PAT system means that two independent process measurements can be taken simultaneously. Optimize and accelerate product development, decrease the time required for QC testing and simplify scale ups while making them more reliable. Gain a deeper knowledge and understanding of the process by observing any aspect of product quality online and in real-time on our intuitive, GAMP-compliant ProcessView control system.

ProcessView control system for easy process control (Precision)

The intuitive touch screen (15" in size) control system with a control panel offers the operator perfect ergonomics. Processes can be controlled easily and reliably.





www.glatt.com

Glatt Systems Pvt. Ltd.

Unit 201-203 & 212-215 (West Wing)
2nd Floor, Nyati Unitree,
Nagar Road, CTS # 1996(B)
Yerwada, Pune -411 006, India

Glatt Systems Pvt. Ltd.

401-404, B Wing,
Dynasty Business Park,
J.B.Nagar,
Andheri-Kurla Road, Andheri (E)
Mumbai -400059, India